

TRUEplus 5-BEVEL STERILE PEN NEEDLES - INSTRUCTIONS FOR USE

INTENDED USE

Sterile, single-use pen needles - TRUEplus 5-BEVEL STERILE PEN NEEDLES are intended for use with a pen injector device for the subcutaneous injection of insulin. **STERILE EO** 

Warning: Remove both the outer protective (clear) cap and the inner protective (color) cap before an injection. If both the outer protective (clear) cap and the inner protective (color) cap are not removed before use, the medication or dose may not be injected, which may result in serious injury or death.

Note: Read the needle gauge and length on the seal to make sure you are using the right pen needle for you. Ask for help if you have impaired vision.


Injection technique: The injection site and technique should be recommended to you by your healthcare professional.



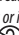


Recommendation: Injection with a skin fold at a 90° angle: It is recommended that you inject into a skin fold vertically (at a 90° angle) to prevent too deep an injection and hitting the muscle tissue. To make a skin fold when injecting, pinch the skin between the thumb and forefinger and do not pull up the muscle. Injection with a skin fold at a 45° angle: This injection technique is only recommended when you have very little subcutaneous fatty tissue. The minimum angle should be 45°. If the angle is less than 45°, the insulin may not reach the subcutaneous fatty tissue and result in wealing or bruising at the injection site. Injection without a skin fold at a 45° angle or a 90° angle: If your layer of subcutaneous fatty tissue is very thick, you can make an injection without a skin fold. When you do not use a skin fold, it is recommended that you inject at a 45° angle and choose a short or medium length needle (6 to 8 mm).



How to get ready for injection: Choose and then disinfect the injection site. Change the injection site often to avoid tissue hardening. You can inject into the fatty tissue of the abdomen, buttocks, or the outer side or front of the upper thigh. Injections into the upper arm or the inner side of the thigh are not recommended, as the risk of injecting insulin into a muscle is high.



DIRECTIONS

- Remove the seal from the bottom of the pen needle.
- Push the pen needle straight on to the pen injector device and twist until it is tight. Do not push the pen needle at an angle.
- Then remove the outer protective (clear) cap and set it aside. Do not discard. You will need the clear cap later to help remove the pen needle from the pen injector device.
- Remove the inner protective (color) cap by pulling it straight out to prevent bending of the needle.
- Prime the pen injector as per the pen's Instructions For Use  to prevent incorrect dosage caused by air in the cartridge or a clogged needle. If no insulin flow is seen after several attempts, replace the outer protective (clear) cap, remove the pen needle from the pen injector device and discard the pen needle. Attach a new pen needle as directed above.
- Set the insulin dose required.
- Inject the selected dose as directed by your healthcare professional; do not bend the needle. Inject slowly (approximately 10 seconds). Do not withdraw the needle immediately after injection. Do not change the direction of the pen needle while it remains in the body, as this can result in bending or breaking of the needle.
- After use, replace the outer protective (clear) cap (only if self-injecting).
- Remove the pen needle from the pen injector device after each injection to prevent air from entering the cartridge and insulin from leaking out. Discard used pen needles according to your local regulations.

Warning: To avoid serious risk of infection, do not use the pen needle after the expiration date  on the packaging has expired. Do not re-use , re sterilize  or share pen needles. Do not use pen needles if the individual packaging is damaged  or if the seal is missing or if the needle appears to be bent or damaged. To avoid accidental needle sticks, never handle someone else's pen needles. Pen needles are single-use only  and used pen needles should be disposed of safely in accordance with your local regulations. Not for use by children without adult supervision.

Caution: Consult your healthcare professional for assistance in choosing the right needle length, injection technique and site. Follow the pen's Instructions For Use . Refer to the Instructions For Use  from the insulin / pen injector manufacturer for specific storage and handling conditions. Seek immediate medical attention if you have any health concerns.

Compatible Injection Pens

INJECTION PEN NAME	MANUFACTURER	INJECTION PEN NAME	MANUFACTURER
apidra SoloStar	Sanofi	NovoLog Mix70/30 FlexPen	Novo
Lantus SoloStar	Sanofi	Victoza Pen	Novo
Humalog KwikPen	Eli Lilly	NovoPen 3	Novo
Humalog Mix 50/50 KwikPen	Eli Lilly	NovoPen 4	Novo
Humalog Mix 75/25 KwikPen	Eli Lilly	NovoPen Echo	Novo
Humulin N Pen	Eli Lilly	NovoPen Junior	Novo
Humulin 70/30 Pen	Eli Lilly	Byetta Pen 5mcg	Amylin
HumaPen Luxura	Eli Lilly	Byetta Pen 10mcg	Amylin
Forteo	Eli Lilly	SymLinPen 60	Amylin
Levemir Flex Pen	Novo	SymLinPen 120	Amylin
Novolog FlexPen	Novo	Autopen	Owen Mumford

Choosing needle length, injection technique and site

AGE	BMI	NEEDLE LENGTH*	
		Without skin fold, straight injection	With skin fold, angled or straight injection
Children < 12 Very Slim Adults	All	4mm, 5mm, 6mm	6mm
Teenagers (12-18 years)	Average	4mm, 5mm, 6mm	8mm
	Overweight (BMI > 25)	4mm, 5mm, 6mm	8mm
Adults	Average	4mm, 5mm, 6mm, 8mm	8mm
	Overweight (BMI > 25)	4mm, 5mm, 6mm, 8mm	8mm, 12mm**, 12.7mm**

*Consult your healthcare professional for assistance in choosing the right needle length, injection technique and injection site.







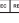
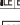

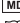


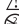




** Recommended for people injecting large dosages or with limited dexterity / strength.

Available TRUEplus 5-BEVEL PEN NEEDLE SIZES

- TRUEplus 5-Bevel pen needles are available in 5 versions, differing in needle diameter and length to meet individual needs.

29G	31G	32G
0.33mm x 12.7mm 29G x 1/2"	0.25mm x 5mm 31G x 3/16"	0.23mm x 4mm 32G x 5/32"
	0.25mm x 6mm 31G x 1/4"	
	0.25mm x 8mm 31G x 5/16"	

SYMBOLS:

 Manufacturer	 Single Use
 Importer	 Lot Number
 Authorized Representative in the European Community	 Sterilized Using Ethylene Oxide
 Use-by Date	 Medical Device
 Caution	 Consult Instructions for Use
 Do not Resterilize	 Do Not Use if Package is Damaged
 Non-Pyrogenic	 Catalogue Number
 Keep Dry	 CE Conformity Mark
	 Keep Away From Sunlight

EC REP **KTR Europe GmbH**
Mergenthalerallee, 77
65760 Eschborn
Germany

Imported by:
Trividia Health UK Limited
27 Old Gloucester Street
London, United Kingdom
WC1N 3AX
UK: 0800 689 5035

CE 1639

TRIVIDIA HEALTH, INC.
2400 N.W. 86TH COURT
FORT LAUDERDALE, FL 33309 U.S.A.
© 12/2024 Trividia Health, Inc. - 56UKT03 Rev. 2